Dama & Sena

SIT BACK AND RELAX





3









Dama

Design, Reliability and Comfort

Treatment beds are an essential part of a practice and should be the focal point of every clinic.

More than just a necessary piece of equipment, Dama represents an investment in quality, technology, and comfort to your patients. It's the opportunity for your clinic to enhance efficiency and expand treatment solutions, differentiating your practice and clearing the way for profit growth.

Design: The design features are augmented by a highly original color concept.

This means that you can explore new avenues in interior design and harmoniously integrate the Dama treatment chair into your clinic. So let yourself be inspired and select your combination.

Reliability: 100% made in Barcelona - Spain. The electronics with SMT technology and the base, seat, back support and leg support are backed up by low tension, silent and reliable equipment.

Comfort: The non-slip upholstery, soft padding and ergonomic lines have been designed to provide maximum comfort to the patient.









Sena

An attractive design with an excellent result

With Dama base the Sena multi-treatment chair offers a new dimension to work: reliability and efficiency at an affordable price.

Easy to operate, safe and economic functioning. This multi-treatment couch allows an easy, versatile work pattern.

Design, comfort and reliability are the main features of this SENA multi-treatment couch. The electronics with SMT (Surface Mounting Technology) and the base, seat, back support and leg support are backed up by low tension, silent and reliable equipment.

The non-slip upholstery, soft padding and ergonomic lines have been designed to provide maximum comfort to the patient.







Details make the difference

Passion for Quality Made in Spain.

It's just a small thing that most people wouldn't notice, but it's a detail worth noting in the quest for outstanding quality and durability.













Cushioned armrests synchronised with the back support, both retractible independently.



Sturdy and comfortable table position to offer a wide range of treatments.



Minimum height from the floor and removable armrests, ensure ease of access for the patient.

Ergonomics and unlimited comfort

The adjustable chair movements, ensure ergonomic working in any position of treatment.

What's your ideal model?

Technical details

Ref.

Ref. with heating

Number of actuators

Power connection

Frecuency

Watt

Standby power

Memories

Warranty

Two-tone colors

Standard color

Movements

Height adjustment

Tilt adjustment

Back adjustment


Legrest adjustment

Link legrest and backrest

Weights

Lifiting capacity


Chair Weight

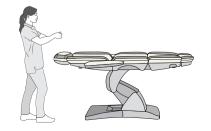


Maximum height of 95 cm allows you to work standing if desired.



Dimensions

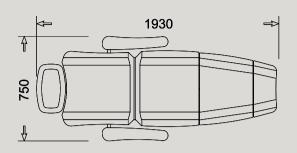

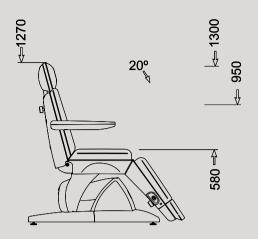

Dama 3 Dama 4 Sena 3 Sena 4

SN-820E	SN-830E	SN-825E	SN-835E
SN-821E	SN-831E		
3	4	3	4
100-240 VAC	100-240 VAC	100-240 VAC	100-240 VAC
50 / 60 Hz			
300 W	350 W	300 W	350 W
Approx. 0.5 W	Approx. 0.5 W	Approx. 0.5 W	Approx. 0.5 W
2	2	2	2
5 years	5 years	5 years	5 years
✓	✓	Х	Х
18	18	white	white

Electrical	Electrical	Electrical	Electrical
Electrical	Electrical	Electrical	Electrical
Electrical	Electrical	Electrical	Electrical
Gas spring	Electrical	Gas spring	Electrical
✓	Memories	✓	Memories

230 Kg	230 Kg	230 Kg	230 Kg
80 Kg.	82 Kg.	82 Kg.	82 Kg.



Patented anti tip system helps to prevent movement when the patient enters or exits the chair.



Powered base, back and tilt for effortless positioning.



Small base footprint which takes up less space in treatment room.



Extendible Head Support

Head support is adjustable both in angle and height, keeping the head steady and a correct position.



Hands free operation

The foot switch is your multifunctional control element for all chair functions. Include home position and two memory positions that you customize to your needs.



Extendible Leg Support

Split extendible legrest, easy to position allows quick and easy adjustments.



Syncro

The legrest raises and lowers in accordance with the backrest. This synchronised movement means the patient is comfortable while the backrest is being adjusted. The Trendelemburg is an asset for safety and comfort.



Massage Headrest

Massage cushion is easily swapped with the standard headrest in just a few seconds. The 3D adjustment reaches a better and more comfortable position.



Thermo heating system

Chair heating system with 3 temperature settings. The chair heating is especially pleasant for the client, and also energy saving, because the temperature in the salon can be set a few degrees lower. (Option not valid for Sena)



Paper Roll Holder

Placed at the rear part of the backrest with quick swap of paper rolls. (max. roll lenght 60 cm, Ø 12,5 cm)



Manicure Bowl

Manicure bowl whith ergonomic support.

Optional Optional Accessoring Figure 1



Protection screen

Protection screen can be attached to either side of the chair.



Casters Set

Locking casters provide mobility without sacrificing valuable floor space. Casters add 3 cm. to the minimum seat height.



Protective Cover

Keeps your chair looking clean and new, protecting the upholstery from the rigors of frequent cleaning.



Magnifying lamp support

Chair side mount magnifying lamp support. Synchronized raising and lowering movement with the chair.



Hand control

The Dama handset has an ergonomic design and size and features large buttons all blending in to a user-friendly handset. Hook for placement on the bed.



Scialytic lamp support

Chair side mount scialytic lamp support. Synchronized raising and lowering movement with the chair.



Wireless Foot Control

Enjoy the freedom of locating the foot control anywhere in the room.

Remove cords and eliminate obstacles.



Stools

Ergonomic for ideal posture.

It's just a small thing that most people wouldn't notice, but it's a detail worth noting in the quest for outstanding quality and durability.



MONZA

Anatomic stool with half moon back rest, upholstered in SKAY, seamless easy to clean with piston controlled height regulation. Soft, double castors mounted on chrome finish, with 5 points of support. Colors can be matched to the range of colors for the treatment beds



COSMOS

Anatomic stool with backrest, upholstered in SKAY, seamless easy to clean with piston controlled height regulation. Soft double wheels mounted on a chrome finish with 5 support points. Colors can be matched to the range of colors for the treatment beds.



CORSA

Ergonomically designed stool to maintain proper posture, upholstered in SKAY seamless, easy to clean with piston controlled height regulation. Soft double wheels mounted on a chrome finish with 5 support points. Colors can be matched to the range of colors for the treatment beds.



CORSA - V

Very healthy, incredibly comfortable "Saddle seat" with backrest. It is possible to adjust the height of the stool through gas spring.



CORSA - L

Very healthy, incredibly comfortable "Saddle seat" with backrest. It is possible to adjust the height of the stool through gas spring.



MONZA - T

It is possible to adjust the height of the stool through gas spring.



Exclusive personality

High-quality skai® material with fine calf leather grain and subdued matte-gloss-effect. The very extensive and trend-orientated color spectrum enables many combinations.

		<u>skai</u>
@	Long-lasting	✓
P	Tear resistant	✓
	Easy to clean	✓
	Abrasion resistance	> 50.000
	Skin irritation and sensitivity ISO 10993-10	No reaction
	Cytotoxicity ISO 10993-5	Complies
	UV resistance	✓
	Composition	90% PVC-compound 10% knitted fabric
*	Made in	Germany

Dama



Sena



The colour samples are indicative and not intended as an exact replica of the actual colour of the product.

The products are subject to modifications and improvements whithout previous notice.



