## PODIATRY





66 When you choose Namrol, you choose something more than a podiatry chair or cabinet.



## MORE THAN TWO DECADES PRODUCING EQUIPMENT OF VERY HIGH QUALITY

Our designers and engineers work together with a common goal, to manufacture the most prestigious and reliable equipment. This has made Namrol one of the leading international companies in the podiatry sector.





















Namrol is the only company prepared to **design, develop, manufacture and market** a complete podiatry clinic.

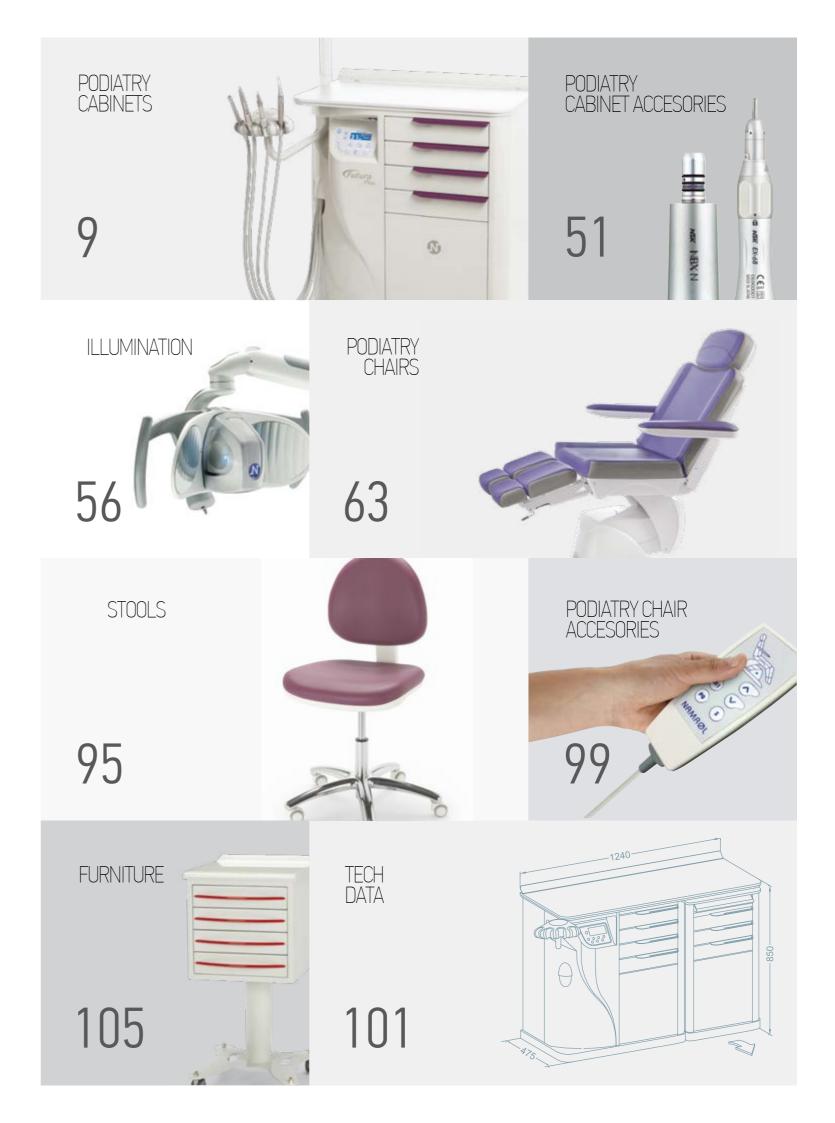





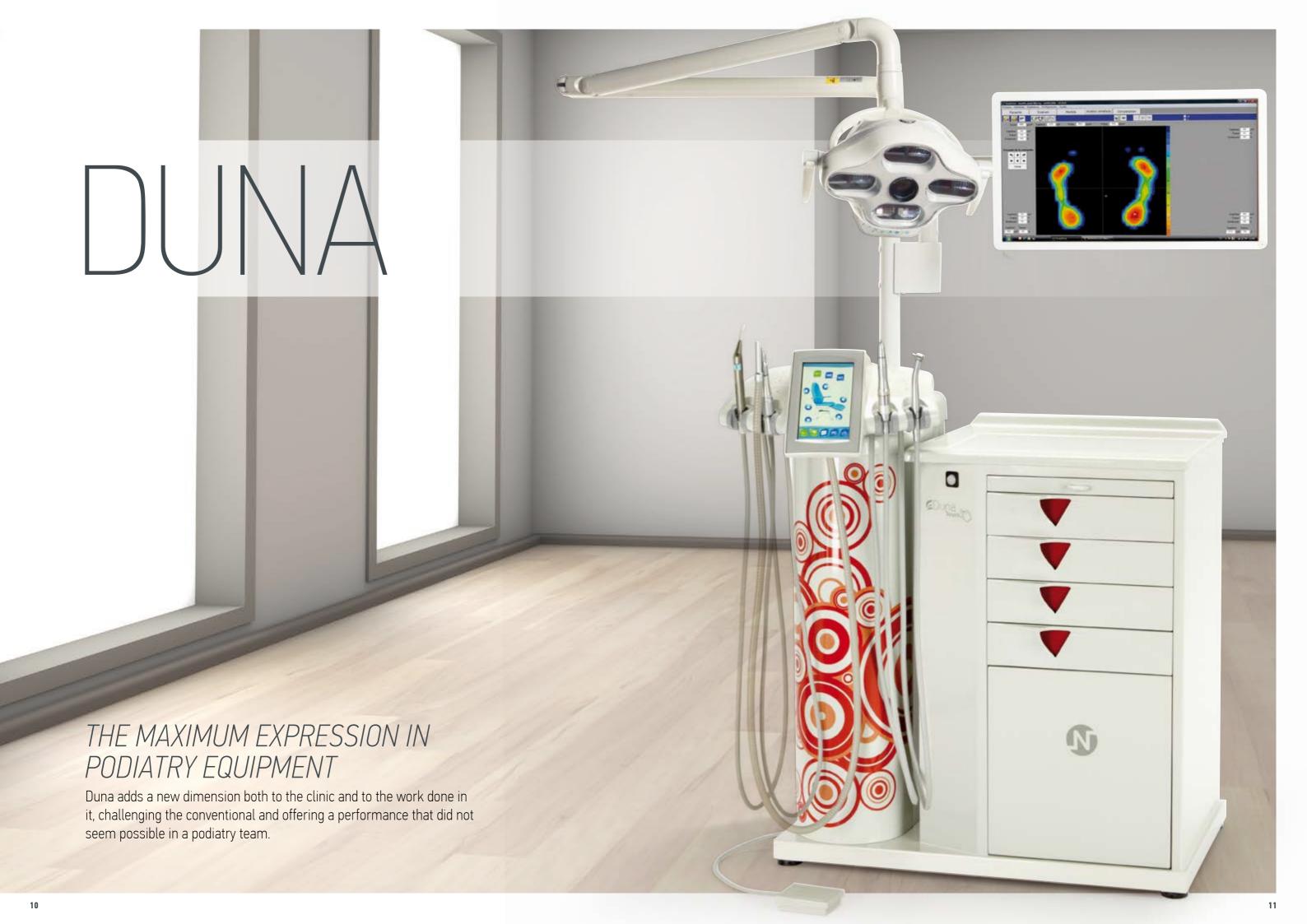



#### Smart or Full Touch console

Duna offers two possibilities, the Smart console with a classic keyboard with LED display and you can also choose the latest generation Full Touch console with 7 "display.

#### Control pedal

All the functions selected in the control panel can be activated by the pedal, thus you will have your hands free. There is also the possibility of having a pedal with integrated speed regulation.

#### Well structured, simple and clear

The display allows the reading and regulation of the revolutions, the direction of rotation of the micromotor and the percentage speed in the case of the turbine. When an instrument is picked up, the settings of the chosen instrument are automatically selected.

#### The best position for each patient

Optionally, an additional module can be incorporated for the wireless control of chair movements from the console.

## INTUITIVE AND EASY-TO-USE CONTROLS

By means of the wide sensitive touch screen you will control your unit even more intuitively, you can visualize and select at any time the parameters of the instruments, the functions of the unit and the chair.

The software warns the user of the need to replace the vacuum bag and other

maintenance processes.



## DRAWER MODULE

Innovative design of handles integrated in the front, fruit of the constant search for functional and ergonomic solutions that guarantee maximum work comfort. Guide system with self-closing compensated. Doors, drawers, front and structure designed to ensure a high level of acoustic comfort. It presents an innovative approach in the choice of solutions and space organization.

#### Soft and silent closing

Namrol incorporates the Softmotion cushioning system. With Sofmotion, the drawers and door close smoothly and quietly, regardless of momentum and weight.


#### Personalization

Generate with colors and shapes an environment where your patients and you feel comfortable. The design indicates to its patients current and progress, key factors to achieve the confidence and fidelity of the patients.

#### Bright and durable surfaces

Materials chosen according to the highest durability requirements in a clinical environment. Resistant while easy to clean.



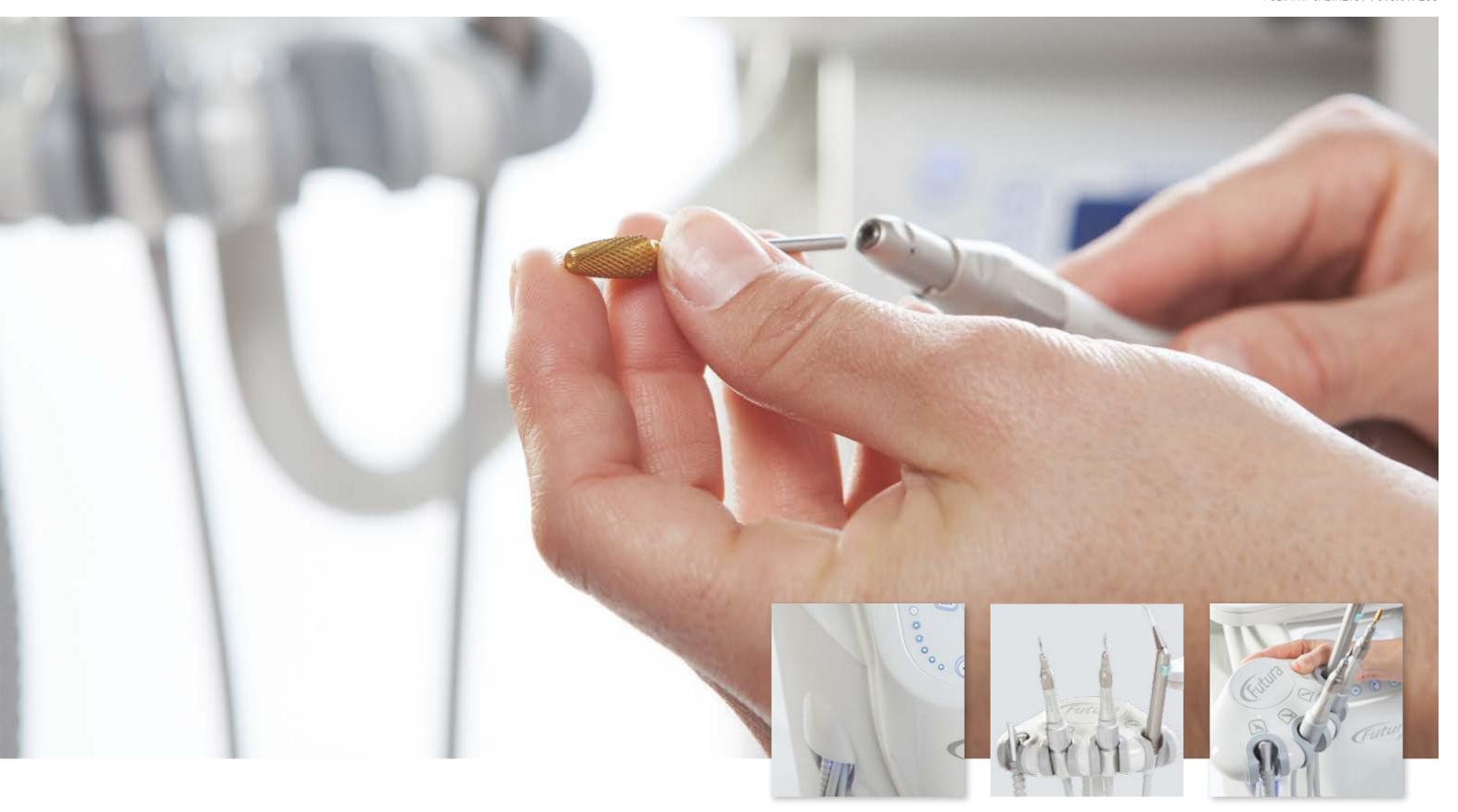

## QUALITY IN ALL DETAILS

The design indicates current and progress, key factors to achieve trust and fidelity of patients. The lacquered pieces of the body of the equipment and the harmonic lines acquire a unique aspect.

#### Bright and durable surfaces

Materials chosen according to the highest durability requirements in a clinical environment. Resistant while easy to clean.

#### Instrument console

The optimal maneuverability of the console of the team favors the maintenance of a correct and natural posture. Stylized and practical, the design of the Futura instrument console allows quick access to instruments. Its configuration of hanging hoses offers a wide range of the instrument making it very manageable.



#### Well structured, simple and clear

The display allows quick reading of the revolutions, the direction of rotation of the micromotor and the percentage speed in the case of the turbine. When an instrument is picked up, the settings of the chosen instrument are automatically selected.

#### **Work Memories**

Frequently used settings (for example: number of revolutions of the instruments) can be stored and retrieved comfortably at any time simply by picking up the instrument.

#### Control pedal

All the functions selected in the control panel can be activated by the pedal, thus you will have your hands free. There is also the possibility of having a pedal with integrated speed regulation.

#### Cleaning the suction tubes

We have designed an innovative automatic cleaning system for the suction pipes that guarantees the efficient operation of the system and a more hygienic treatment environment. The circuit is automatically cleaned every time we hang the instrument.

## INTUITIVE CONTROL PANEL

The control keyboard manages all the functions of the unit, achieving maximum simplicity and total information on the LCD screen. The keys ensure a simple sanitization at the same time as a quick adjustment of the instruments and the functions of the unit with automatic selection when the instrument is picked up. Easy pre-selection reading of the set speed of the micromotor in rpm. and the direction of rotation. It incorporates two work memories for each instrument.



#### Solid surface

The table top manufactured in non-porous Quirofibra helps to avoid the penetration of liquids and facilitates cleaning and disinfection.

#### Automatic opening without handle

The door without handle is opened without effort by a slight pressure on the front. To close the door, simply press it lightly.

#### Easy access

Access to the dust deposit is made from the front door to facilitate the replacement of the bag. The bottle of liquids incorporates an automatic system of depressurization to facilitate the filling.

## SOFT AND SILENT CLOSING

Namrol incorporates the Softmotion cushioning system. With Sofmotion, the drawers and door close smoothly and quietly, regardless of momentum and weight. Door, drawers, front and structure designed to ensure a high level of acoustic comfort. It presents an innovative approach in the choice of solutions and space organization.







## QUALITY IN ALL DETAILS

The design conveys current and progress, key factors to achieve the trust and fidelity of patients. The lacquered parts of the body of the equipment and the harmonic lines acquire a unique aspect..

#### Left-handed version

Same features and features as the right-handed version

#### Intuitive control panel

The control keyboard manages all the functions of the unit, achieving maximum simplicity. The keys ensure a simple sanitization at the same time as a quick adjustment of the instruments and the functions of the unit with automatic selection when the instrument is picked up. It incorporates two work memories for each instrument.

#### Instrument console

The optimal maneuverability of the console of the team favors the maintenance of a correct and natural posture. Stylized and practical, the design of Stylo's instrument console allows quick access to instruments. Its configuration of hanging hoses offers a wide range of the instrument making it very manageable.

#### Cleaning the suction tubes

We have designed an innovative automatic cleaning system for the suction pipes that guarantees the efficient operation of the system and a more hygienic treatment environment. The circuit is automatically cleaned every time we hang the instrument.

29



#### Solid surface

The table top manufactured in non-porous Quirofibra helps to avoid the penetration of liquids and facilitates cleaning and disinfection.

#### Automatic opening without handle

The door without handle is opened without effort by a slight pressure on the front. To close the door, simply press it lightly.

#### Easy access

Access to the dust deposit is made from the front door to facilitate the replacement of the bag. The bottle of liquids incorporates an automatic system of depressurization to facilitate the filling.

## SOFT AND SILENT CLOSING

Namrol incorporates the Softmotion cushioning system. With Sofmotion, the drawers and door close smoothly and quietly, regardless of momentum and weight. Door, drawers, front and structure designed to ensure a high level of acoustic comfort. It presents an innovative approach in the choice of solutions and space organization.





## DRAWER MODULE

Innovative design of handles integrated in the front, fruit of the constant search for functional and ergonomic solutions that guarantee maximum work comfort. Guide system with self-closing compensated. Doors, drawers, front and structure designed to ensure a high level of acoustic comfort. It presents an innovative approach in the choice of solutions and space organization.

#### Compact volumes

Maximum comfort of work. Guide system with self-closing compensated. Doors, drawers, front and structure designed to ensure a high level of acoustic comfort. It presents an innovative approach in the choice of solutions and space organization.

#### Few keys, many features

Namrol has adopted the technology of simplicity, which allows many features and information to be obtained with few keys. The keys under the membrane allow to manage all the main functions, always keeping the state of the instruments under control and in a simple way.

## Instruments that can be handled almost effortlessly

The adjustable instrument console has been conceived for the professional who works with superior recovery instruments. The rotation of the support and the movement of the hoses on guides, allows positioning the console in the right place.





### FUNCTIONALITY AND VERSATILITY

Under the simplicity of its design, it hides a treatment unit of great technical complexity, achieved by its monoblock metal construction and the best ergonomic solutions, to ensure the podiatrist high performance and optimal functionality in your clinic.

The optimal maneuverability of the console of the team favors the maintenance of a correct and natural posture.

#### Well structured, simple and clear

The control keyboard manages all the functions of the unit. The keys ensure a simple sanitization at the same time as a quick adjustment of the instruments and the functions of the unit with automatic selection when the instrument is picked up. It incorporates two work memories for each instrument.

#### Solid surface

The table top manufactured in non-porous Quirofibra helps to avoid the penetration of liquids and facilitates cleaning and disinfection. It incorporates a handle for the movement of the unit.

#### Mobility

Optionally Activ can incorporate a system of wheels with individual locking, facilitating movement in several rooms of the clinic.

#### Connections

On the back it incorporates a 230V connection, ideal for the incorporation of magnifying lamps.





#### Smart management concept

The membrane keyboard technology allows intuitive and efficient instrument handling. The Epsylon Plus keyboard also incorporates the chair controls.

#### Bright and durable surfaces

Materials chosen according to the highest durability requirements in a clinical environment. Resistant while easy to clean.

#### Auxiliary tray

The removable auxiliary tray, located in the lower part of the console, allows to have the instruments near the work area, saving unnecessary movements.

#### Easy and logical

The auxiliary tray incorporates a glove dispenser at the bottom.

## FLEXIBILITY AND ADAPTABILITY

Thanks to the concept of arm mounted on the chair, the podiatrist is in optimal conditions to work in a small space. The strong point of Epsylon Plus is its mobility, flexibility and adaptability to all operational situations.



#### External deposit

Access to the independent water supply tank is carried out from the outside to facilitate its replacement or simple filling.

#### Closeness to the patient

In all situations, comfort has a constant and shared role.

#### Compact and balanced in size

The instrument console can be placed near the patient's foot, which facilitates the reach of instruments and the keyboard.

#### Left-handed version

Same features and features as the right handed version.

By means of the wide sensitive touch screen you will control your unit even more intuitively, you can visualize and select at any time the parameters of the instruments, the functions of the unit and the chair.

The software warns the user of the need to replace the vacuum bag and other

maintenance processes.





#### Adjustable instrument console

Console with double articulation able to rotate around its axis to reach the most ideal position in relation to the type of treatment performed, providing maximum freedom of movement. The functional arrangement of the hanging hose instruments in the console ensures an immediate and comfortable grip. The successive return of them occurs intuitively and precisely, maintaining a quick and easy work.

#### Smart management concept

The membrane keyboard technology allows an intuitive and efficient operation. When an instrument is picked up, the settings of the chosen instrument are automatically selected. The keyboard layout allows a quick adjustment of the instruments and functions of the unit with reading of the pre-selection of the set speed of the micromotor in rpm. and direction of rotation.

## Auxiliary tray for an agile work in its daily use

The removable auxiliary tray, located on the side of the console, allows to have the instruments near the work area, saving unnecessary movements.

#### External deposit

Access to the independent water supply tank is carried out from the outside to facilitate its replacement or simple filling.

## FLEXIBILITY AND ADAPTABILITY

Thanks to the concept of arm mounted on the chair, the podiatrist is in optimal conditions to work in a small space. The strong point of Epsylon is its mobility, flexibility and adaptability to all operational situations.





#### **GT40 Micromotor**

The GT40 electric micromotors offer high performance. One of the most effective and light micromotors in podiatry. Work without vibrations up to 40,000 rpm., With a high torque at low revolutions.

#### MC2IR Micromotor

Electric micromotor with internal spray, reliable and resistant. More than remarkable performances thanks to a speed of 40,000 rpm.

Electric micromotor with internal spray, reliable and resistant. Titanium outer body, ultra-compact and lightweight. Minimum generation of heat. More than remarkable performances thanks to a speed of 40,000

**NBX Micromotor** 

#### MCX Micromotor

Brushless micromotor (brushless), maintenance-free with a large torque at low and high speeds between 1000 and 40,000 rpm. Thanks to its reduced size by 30%, it adapts perfectly to the hand and thus guarantees an optimal balance with a perfect weight distribution.

#### PM70 by NSK Hand piece

The PM70 handpiece is specially designed for use in podiatry. The sophisticated technology of ceramic bearings, increases the useful life and the good operation. No vibrations and very silent.

### X65 Hand piece

Straight hand piece with internal spray, titanium body with Durocoat coating. No vibrations and very silent. Optionally you can incorporate LED lighting.

#### PM60 by Bien-Air Hand piece

Straight hand piece with internal spray, designed to withstand the great tensions in repeated lateral supports and maintain an incomparable precision. No vibrations and very silent.

#### Suction cannula

Made of anodized aluminum, it provides lightness and strength. For a clean and dust-free job.

#### Black Pearl turbine

Unsurpassed quality / price ratio. Designed to offer a power 310,000 rpm. and a robustness more than remarkable. As for the change of milling cutter, it is as simple thanks to the push button integrated in the head. Direct connection to 4-hole hose.

#### Spraytec 3F syringe

3 functions stainless steel syringe. Guaranteed operation thanks to the push button system that separates the drive from the air and water.



#### Led lamp with integrated HD camera

LED lamp with integrated video camera with FullHD resolution (1920x1080px) with autofocus, allowing you to see the work area with a 30x optical zoom.

The potential of this tool, applied to the field of podiatry are remarkable: see the detailed image of the intervention, highlighting the details barely visible to the naked eye, increases the visual capabilities of the professional and relieves eye fatigue, leading to an important advance in the methodology of work. The use of the camera facilitates communication with the patient. In the diagnosis or treatment phase, the pathologies can be clearly shown, which helps an immediate understanding of the problem. During the treatment phase, the podiatrist can follow in real time or record the evolution of the intervention, after which it will be possible to evaluate the results of the work carried out.

#### 19 "or 22" monitor with articulated arm

By choosing the double articulation arm, the professional can move the monitor and orient the screen within the operative area according to their own needs, improve communication with the patient, show the treatments performed and document the clinical history with images. The monitor is predisposed for wired connection to a PC, via HDMI connector.

#### Integrated removable keyboard

Provided with USB connection; It allows access to all the functions of the cabinet software without leaving the patient. It is resistant and is hidden in the unit, thus leaving a tabletop for free work.

## MULTIMEDIA FOR YOUR CLINIC

The multimedia accessories for Namrol units, open a new vision that meets all expectations regarding modern communication for consultation. They help reduce diagnostic times, show and share the images with the patient, convince them to carry out the treatment

54 5!

# ALYA PODO











3.000 a 40.000 Lux

5000 K (day light)

LED

23VA

#### Hygiene

Light intensity

Illumination

Absorbed power

Color temperature

The two handles are removable and can be sterilized.

**TECHNICAL CHARACTERISTICS** 



#### 3D rotation

The perfect balance, combined with the exclusive 3-D rotation system of the head on the optical axis, ensures the best orientation of the light source.



#### Perfect adjustment

The new electronic control of the illumination, activated with joystick, regulates the light in a continuous way, keeping the temperature of the color constant.



From the great experience of Namrol in podiatry Alya Podo was born, a lamp with a larger illuminated field of 30x20 cm. that make it the ideal lamp for podiatry.

Thanks to the exclusive optical projector of reflected light, through the 258 small reflecting cells of each of the two parabolas, the direct vision of the 2 leds is avoided

Long lifetime and reduced power consumption thanks to advanced diode technology with fanless cooling.



#### Hygiene

The two handles are removable and can be sterilized.



#### Projection

The precise definition of the light source eliminates the risk of glare.



#### Perfect adjustment

The new electronic control of the illumination, activated with joystick, regulates the light in a continuous way, keeping the temperature of the color constant.

## TECHNICAL CHARACTERISTICS Polyv

Light intensity	35.000 Lux
Color temperature	5000 K (day light)
Illumination	LED
Absorbed power	9VA

### Polyvalent and functional

The new lamp, for those looking for high performance LED at a competitive price. Thanks to the exclusive optical projector of reflected light, through the 258 small reflecting cells of each of the two parabolas, the direct vision of the 2 LEDs is avoided, creating a totally uniform illuminated field of 175 x 100 mm.

Long lifetime and reduced power thanks to advanced diode technology with fanless cooling.





The halogen light source is diffused through a double rear parabolic glass that guarantees a light without shadows.




Functionality and flexibility

Articulated arms that allow to easily orient the light, with smooth and precise movements.



**Control electronics** 

On and off with switch located in the head, with forced ventilation system regulated by microprocessor.

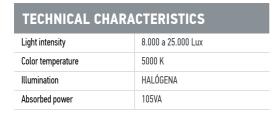


#### Very versatile

EDI is an advanced design lamp, which guarantees excellent illumination of the work area. The design combines ergonomics and simplicity of intervention, with smooth and rounded surfaces to facilitate cleaning.



Intelligent control

On and off synchronized with the head, or by the push button integrated in the arm. Automatic storage of lighting.

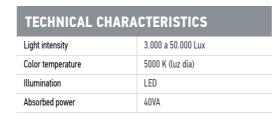



Light support

ALYA + THEIATECH

Integrated in the back arm of Alya, THEIATECH has been conceived to work together with the operating light.



LED technology

La fuente luminosa de THEIATECH está formada por 54 LEDS que, a través de un específico difusor, producen un cono de luz en el entorno circundante.

#### The artificial light closest to the sun

THEIATECH has been conceived to create a luminous zone of transition between the work area and the instrumentation area, balancing the illumination of the areas of the podiatric cabinet and reducing visual fatigue. Thanks to a complete illumination of the working area of between 15,000 and 50,000 lux, a constant focus of the eye is avoided and the iris contracts and dilates repeatedly. This overload of the visual apparatus is detrimental to the health of the eyes and can impair good health at work.

# AEROBOX

POWER AND SILENCE UNITED



ONE OF THE QUIETEST COMPRESSORS ON THE MARKET





#### Powerful and durable

Turbine suction motor, with maximum dust extraction efficiency.

TECH. CHARACTERISTICS

Pressure

Noise level

Absorbed power

Flow

1300 mm H20

1.350 l/min.

56 dB.

3,1A



#### Adaptable

Optionally, Aerobox can be equipped with a cannula aspiration system.



#### Clean

The exchangeable bag system provides cleaning and convinces by a hygienic filter system.



Let yourself be convinced by a more powerful and at the same time more silent aspiration. Due to its small size, it is ideal for new installations or the update of existing installations. The use of a brushless suction motor, cancels mechanical friction and provides significantly longer engine operation. Compact steel structure with foamed interior, soundproof and antibacterial.



#### Hygiene

Quiro-Air Plus incorporates as standard a submicron filter (0.01 µm) filtering any particle in the supplied air.



#### Autonomous

NAMRØL

•

Due to its low level of sound, it can be operated autonomously thanks to its simple plug-and-play features.



One of the main features is space saving due to its compact design.

NAMROL

TECH. CHARACT.	QUIRO-AIR	QUIRO-AIR PLUS
Air flow aspirated	65 l/min.	100 l/min.
Noise level	42 dB.	42 dB.
Absorbed power	0,50 kw.	0,75 kw.
Filtering	5 μm	0,01 µm

#### Quiro-air silent compressor

Namrol's new silent compressor (patent U9701105) produces quality compressed air, ideal to meet the demands of the podiatric clinic.

It has a direct drive electric motor with built-in overload protection for maintenance-free operation

Manufactured with high quality components and materials, it is a reliable source of high quality air. These compressors are compact and flexible, and are available in two versions depending on size.







## FUNCTION AND SHAPE MEETS TOGETHER

A design with tenuous vintage touches that is accompanied by a design that puts ergonomics and comfort beforehand to enhance the welcoming feeling and well-being of the patient.

#### Smart management concept

The control pedal has buttons for direct action to the automatic memories and zeroing, which together with the easy-memory system for quick memory and recovery of work positions, maximizes functionality and simplicity of work.

#### Two-piece double leg loops

Leg with double support area for heel and calf. Extensible, with lateral opening, rise and fall compensated by gas piston.

#### Integrated headrest

You can easily move vertically to adapt it to people of different heights.

### Color combination

The attractive contrast upholstery gives Nova padding a first class look.

## SMART AND ADAPTABLE BACK

The seat is the starting point from which the patient rests on a completely ergonomic backrest, which facilitates total relaxation and comfort during the



### CHARACTERISTICS

TREATMENT CHAIR	NOVA 3	NOVA 4
Motors	3	4
Weight	96 Kg.	97 Kg.
Maximum lifting capacity	220 Kg.	220 Kg.
Upholstered thickness	7 - 10 cm.	7 - 10 cm.
Available colours	50 Combinations	50 Combinations
High gloss ABS housings	✓	✓
Turn 360° with lock	✓	✓
Extendable leg loops with side opening	✓	✓
2 Automatic memories	✓	✓
Automatic zeroing	✓	✓
Height movement	Electric	Electric
Movement of the seat - trendelemburg	Electric	Electric
Movement of the backrest	Electric	Electric
Movement of the leg	Gas shock absorber	Electric
Ergonomic head	Integrated	Integrated
Removable rigid tray	Optional	Optional
Arm with lamp holder	Optional	Optional



LINAK む

Shock absorbers

skai

**STABILUS** 





Swivel base The base allows a 360o rotation for optimal work

positioning.



The best position Multiple movements to get the perfect position for each treatment.



Foot treatment With a position of maximum height of 1,300

mm. You can also work standing up.



Adaptable leg loops Leg extensible, folding and with lateral opening.



Anti-roll system Patented anti-roll system to prevent movements during patient entry or exit.



Space for the legs The reduced design of the base provides more free space for the legs.

Upholstered arms with lowered synchronization with the movement of the backrest, folding independently.

The swivel base, the minimum height to the floor and the folding arms, facilitate a comfortable access to the patient.

Anatomical backrest with low hardness foam, for a comfortable and relaxed treatment of patients.





### Smart management concept

The control pedal has buttons for direct action to the automatic memories and zeroing, which together with the easy-memory system for quick memory and recovery of work positions, maximizes functionality and simplicity of work.

### Premium padding

The padding with non-deformable polyurethane foam convinces not only for its high level of comfort, but also for its durability and good appearance over time.

### Cushioned leg

Single leg, extendable with rise and fall compensated by double gas piston.

## MAKE YOUR PATIENTS FEEL COMFORTABLE

Perfect combination of form and function. Its ergonomic design conforms to the shape of the patient's body and provides a firm and comfortable support. The seamfree upholstery provides a pleasant feel to the touch and can be easily cleaned and disinfected.

# INDIVIDUAL **POSITIONING** AND ACCESS COMFORTABLE OF THE PATIENT

The unique collapsible leg facilitates patient access and provides a comfortable vertical sitting position for consultation exams and after a treatment.



### **CHARACTERISTICS**

TREATMENT CHAIR	SIGMA 3	SIGMA 4
Motors	3	4
Weight	95 Kg.	97 Kg.
Maximum lifting capacity	220 Kg.	220 Kg.
Upholstered thickness	7 cm.	7 cm.
Available colours	18	18
High gloss ABS housings	✓	✓
Turn 360o with lock	✓	✓
Extendable leg loops with side opening	✓	✓
2 Automatic memories	✓	✓
Automatic zeroing	✓	✓
Height movement	Electric	Electric
Movement of the seat - trendelemburg	Electric	Electric
Movement of the backrest	Electric	Electric
Movement of the leg	Gas shock absorber	Electric
Ergonomic head	Optional	Optional
Removable rigid tray	Optional	Optional
Arm with lamp holder	Optional	Optional

LINAK む

Shock absorbers



**STABILUS** 





Upholstered arms with lowered synchronization The swivel base, the minimum height to with the movement of the backrest, folding the floor and the folding arms, facilitate a comfortable access to the patient.

Anatomical backrest for a comfortable and relaxed treatment of patients. Adaptation of the chair to the height of the patient, which allows an individual and very comfortable positioning of the patient.



Swivel base The base allows a 360o rotation for optimal work positioning.



The best position Multiple movements to get the perfect position for each treatment.



Foot treatment With a position of maximum height of 1,300 mm. You can also work standing up.



Adaptable leg Extendable leg and folding.



Anti-roll system Patented anti-roll system to prevent movements during patient entry or exit.



Space for the legs The reduced design of the base provides more free space for the legs.

independently.





### Smart management concept

The control pedal has buttons for direct action to the automatic memories and zeroing, which together with the easy-memory system for quick memory and recovery of work positions, maximizes functionality and simplicity of work.

### Premium padding

The padding with non-deformable polyurethane foam convinces not only for its high level of comfort, but also for its durability and good appearance over time.

### Cushioned leg loops

Double leg, extendable with rise and fall compensated by double gas piston.

### FIRST CLASS TREATMENTS

Omega stands out for the perfect combination of form and function. Its ergonomic design conforms to the shape of the patient's body and provides a firm and comfortable support. The seam-free upholstery provides a pleasant feel to the touch and can be easily cleaned and disinfected.

# LUMBAR SUPPORT **ASYMMETRIC**

The shapes have been conceived from the study of ergonomics, physiognomy and kinematics of the human body for optimal posture during treatment.



TREATMENT CHAIR	OMEGA 3	OMEGA 4
Motors	3	4
Weight	95 Kg.	97 Kg.
Maximum lifting capacity	220 Kg.	220 Kg.
Upholstered thickness	7 cm.	7 cm.
Available colours	18	18
High gloss ABS housings	✓	✓
Turn 360° with lock	✓	✓
Extendable leg loops with side opening	✓	✓
2 Automatic memories	✓	✓
Automatic zeroing	✓	✓
Height movement	Electric	Electric
Movement of the seat - trendelemburg	Electric	Electric
Movement of the backrest	Electric	Electric
Movement of the leg	Gas shock absorber	Electric
Ergonomic head	Optional	Optional
Tray with removable leg protector	Optional	Optional
Removable rigid tray	Optional	Optional
Arm with lamp holder	Optional	Optional

LINAK<sup>®</sup>む

Shock absorbers











Swivel base The base allows a 360o rotation for optimal work positioning.



The best position Multiple movements to get the perfect position for each treatment.



Foot treatment With a position of maximum height of 1,300 mm. You can also work standing up.



comfortable access to the patient.

Anatomical backrest for a comfortable and relaxed treatment of the patient.



Adaptable leg loops Leg extensible, folding and with lateral opening.



Anti-roll system Patented anti-roll system to prevent movements during patient entry or exit.



Space for the legs The reduced design of the base provides more free space for the legs.





## RELIABILITY AND COMFORT

Aesthetics with faint vintage touches that is accompanied by a design that puts ergonomics and comfort before the patient's welcoming sensation and comfort.

### Smart management concept

The control pedal has buttons for direct action to the automatic memories and zeroing, which together with the easy-memory system for quick memory and recovery of work positions, maximizes functionality and simplicity of work.

### Two-piece double leg loops

Leg with double support area for heel and calf. Extensible, with lateral opening, rise and fall compensated by gas piston.

### Integrated headrest

You can easily move vertically to adapt it to people of different heights.

### Color combination

The attractive contrast upholstery gives the Elyse padding a first class look.

# SMART AND ADAPTABLE BACK

The seat is the starting point from which the patient rests on a completely ergonomic backrest, which facilitates total relaxation and comfort during the



TREATMENT CHAIR	ELYSE 3	ELYSE 4
Motors	3	4
Weight	96 Kg.	97 Kg.
Maximum lifting capacity	220 Kg.	220 Kg.
Upholstered thickness	7 - 10 cm.	7 - 10 cm.
Available colours	50 Combinations	50 Combinations
High gloss ABS housings	✓	✓
Extendable leg loops with side opening	✓	✓
2 Automatic memories	✓	✓
Automatic zeroing	✓	✓
Height movement	Electric	Electric
Movement of the seat - trendelemburg	Electric	Electric
Movement of the backrest	Electric	Electric
Movement of the leg	Gas shock absorber	Electric
Ergonomic head	Integrated	Integrated
Removable rigid tray	Optional	Optional
Arm with lamp holder	Optional	Optional

TOP QUALITY COMPONENTS

Motors and control

Shock absorbers











Upholstered arms with lowered synchronization with the movement of the backrest, folding independently.

to the patient.

The minimum height to the ground and the folding arms, facilitate a comfortable access

Anatomical backrest with low hardness foam, for comfortable and relaxed treatment of patients. Anatomical backrest with low hardness foam, for a



Stability The broad base allows all

types of treatments with complete stability.



The best position Multiple movements to get the perfect position

for each treatment.



Foot treatment

With a position of maximum height of 1,300 mm. You can also work standing up.



Adaptable leg loops Leg extensible, folding and with lateral opening.



Anti-roll system Patented anti-roll system to prevent movements during patient entry or exit



Space for the legs The reduced design of the base provides more free space for the legs.





## RELIABILITY AND COMFORT

Aesthetics with faint vintage touches that is accompanied by a design that puts ergonomics and comfort before the patient's welcoming sensation and comfort.

### Smart management concept

The control pedal has buttons for direct action to the automatic memories and zeroing, which together with the easy-memory system for quick memory and recovery of work positions, maximizes functionality and simplicity of work.

### Two-piece double leg loops

Leg with double support area for heel and calf. Extensible, with lateral opening, rise and fall compensated by gas piston.

### Integrated headrest

You can easily move vertically to adapt it to people of different heights.

### Color combination

The attractive contrast upholstery gives the Elyse padding a first class look.





# FAST MOVEMENTS, SILENT AND RELIABLE

An armchair, for professional use, that incorporates aesthetics and versatility. Penta is committed to design, soft and rounded shapes that suggest comfort and make it perfect for any current clinic.



### **CHARACTERISTICS**

TREATMENT CHAIR	PENTA
Motors	1
Weight	85 Kg.
Maximum lifting capacity	220 Kg.
Upholstered thickness	7 cm.
Available colours	16
High gloss ABS housings	√
Turn 360° with lock	√
Extendable leg loops with side opening	√
Automatic zeroing	√
Height movement	Electric
Movement of the seat - trendelemburg	Compensated with the rise
Movement of the backrest	Gas shock absorber
Movement of the leg	Gas shock absorber
Ergonomic head	Optional
Tray with removable leg protector	Optional
Removable rigid tray	Optional
Single leg	Optional

Motors and control

Shock absorbers











The swivel base, the minimum height to the floor and the folding arms, facilitate a comfortable access to the patient.

Anatomical backrest for a comfortable and relaxed treatment of patients. Adaptation of the chair to the height of the patient, which allows an individual and very comfortable positioning of the patient.



Swivel base The base allows a 140o turn for optimal work positioning.



Adaptable leg loops Leg extensible, folding and



The best position

Multiple movements to get the perfect position for each treatment.



Foot treatment

With a position of maximum height of 1,100 mm. You can also work standing up.



Anti-roll system Patented anti-roll system to prevent movements during patient entry or exit



**Space for the legs**The reduced design of the base provides more free space for the legs.

Ergonomic folding arms

independently.





### Cosmos

The ergonomics of the Cosmos chair has been designed to meet the needs of the clinical user.

Therefore, among its features are the permanent contact backing, skay upholstery, without easy cleaning seams and height adjustment by gas piston. Soft, double wheels, mounted on a chromed base with 5 support points.

### Monza

Monza responds to the most demanding clinics since it brings together everything a cabinet might need. Half-moon backrest, skay upholstery, easy cleaning and adjustment of the height by gas piston. Soft, double wheels, mounted on a chromed base with 5 support points.

### Corsa

Ergonomic stool designed to keep the spine straight, upholstery in skay, without easy cleaning seams and height adjustment by gas piston. Soft, double wheels, mounted on a chromed base with 5 support points. The color can be combined with the range of colors of armchairs.



# ACCESORIOS SILLONES Complemente su sillón con los accesorios adicionales



### Controller

Manual control system with easy-to-clean keypad. (only valid for armchairs of 3 and 4 engines)



Roll holder

Ergonomic support for easy roll placement (maximum length 59 cm.)



Lamp holder magnifying glass

Ergonomic support for easy roll placement (maximum length 59 cm.)



Support for light bulb

System attached to the chair for mounting lamps with articulated arm.



Tray with protector

Waste collection trays with leg protector and quick fastening system.



Tray

Waste collection trays with leg protector and quick fastening system.



**Ergonomic head** 

Maximum comfort and support for your patients. Adjustable in height






## AUXILIARY MODULES

Linear forms in combination with the curves and the design of the countertop offer an innovative line, and refined style, manufactured in Quirofibra® non-porous and without joints, everything is designed for easy cleaning and maximum germ reduction.






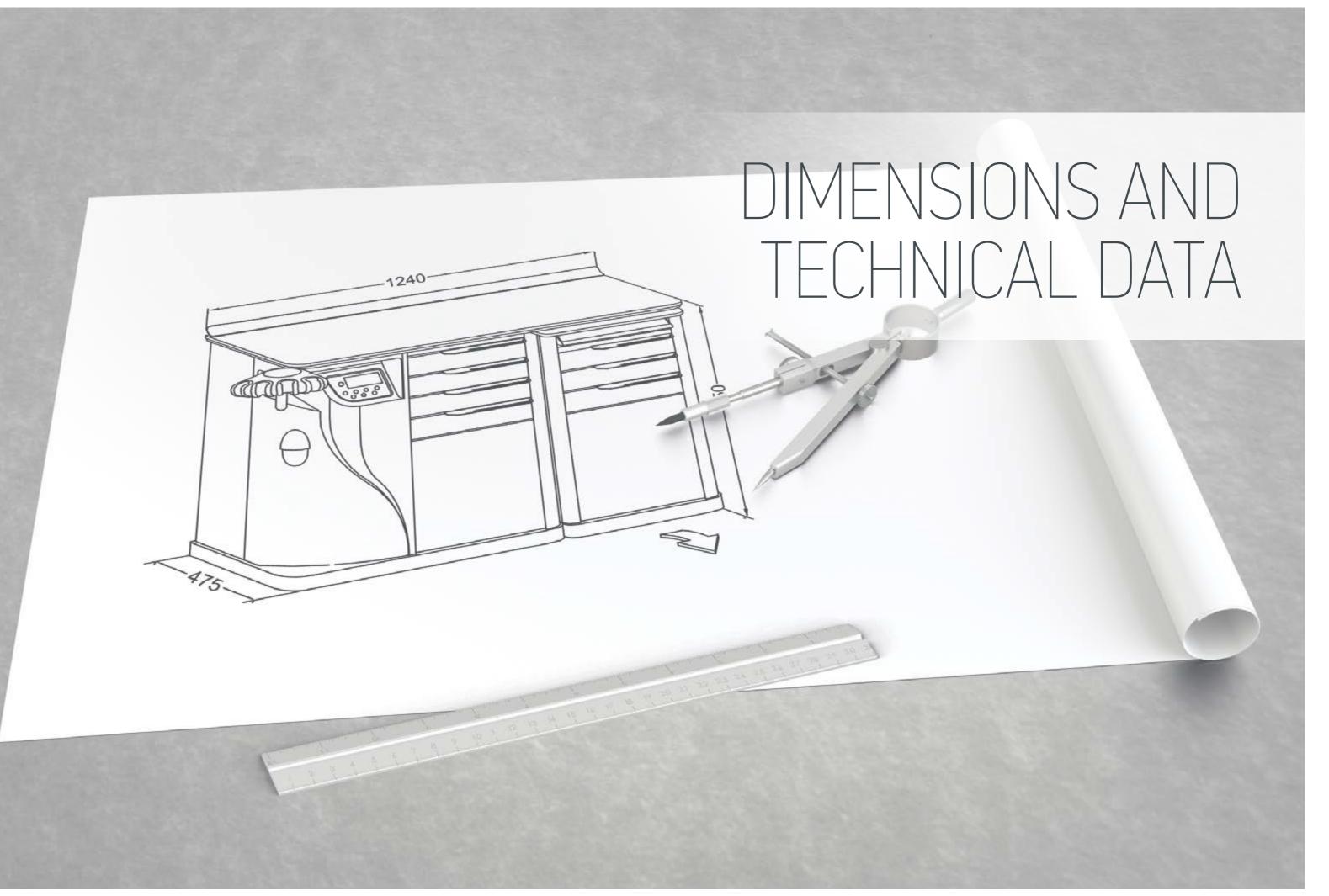

## DUNA

**POWER SUPPLY:** 230 V. 50/60 HZ ABSORBED POWER: 4 A.

MAXIMUM NUMBER OF INSTRUMENTS: 4
INTERNAL PRESSURE: AIR 3,5 BAR / 5,5 BAR. WATER SPRAY 2,5 BAR

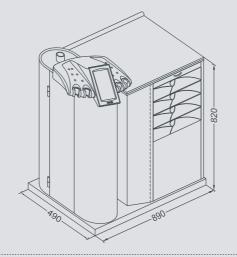


## FUTURA PLUS

POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 4 A.

WEIGHT: 85 KG.

MAXIMUM NUMBER OF INSTRUMENTS: 4
INTERNAL PRESSURE: AIR 3,5 BAR / 5,5 BAR. WATER SPRAY 2,5 BAR

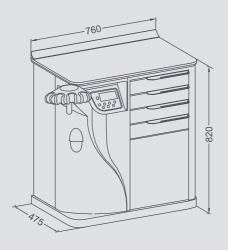


## **FUTURA XL**

**ALIMENTACIÓN ELÉCTRICA:** 230 V. 50/60 HZ **ABSORBED POWER:** 4 A.

WEIGHT: 105 KG.

MAXIMUM NUMBER OF INSTRUMENTS: 4 INTERNAL PRESSURE: AIR 3,5 BAR / 5,5 BAR. WATER SPRAY 2,5 BAR

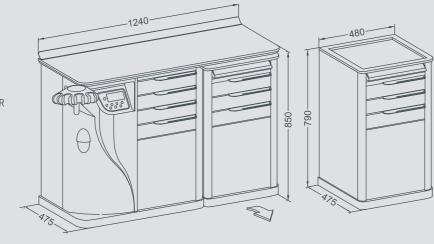


## **STYLO**

POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 4 A. WEIGHT: 85 KG.

MAXIMUM NUMBER OF INSTRUMENTS: 3


INTERNAL PRESSURE: AIR 3,5 BAR / 5,5 BAR. WATER SPRAY 2,5 BAR

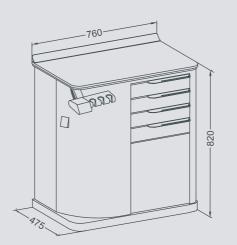

## MAXIMA


POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 4 A.

MAXIMUM NUMBER OF INSTRUMENTS: 5
INTERNAL PRESSURE: AIR 3,5 BAR / 5,5 BAR. WATER SPRAY 2,5 BAR

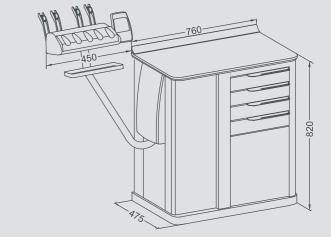

## **ACTIV**

POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 3 A. WEIGHT: 57 KG. MAXIMUM NUMBER OF INSTRUMENTS: 3
INTERNAL PRESSURE: AIR 3 BAR. WATER SPRAY 2,5 BAR

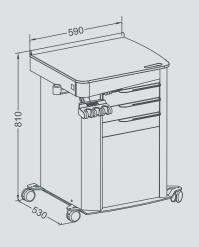


## **EPSYLON PLUS**


(Adaptable to Omega and Sigma)

POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 4 A.

MAXIMUM NUMBER OF INSTRUMENTS: 4

INTERNAL PRESSURE: AIR 3,5 BAR / 5,5 BAR. WATER SPRAY 2,5 BAR

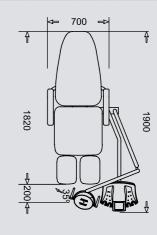

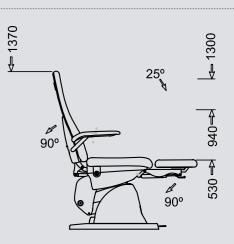

## **EPSYLON**


(Adaptable to Beta and Delta)


POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 4 A.

WEIGHT: 90 KG.

MAXIMUM NUMBER OF INSTRUMENTS: 4


INTERNAL PRESSURE: AIR 3,5 BAR / 5,5 BAR. WATER SPRAY 2,5 BAR

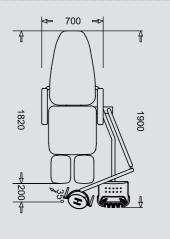


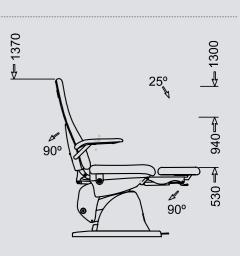

## NOVA


POWER SUPPLY: 230 V. 50/60 HZ
ABSORBED POWER: 3 A.
WEIGHT: 90 KG.
MOVEMENT: ELECTROMECHANICAL
MOTORS AND ELECTRONICS: LINAK
MAXIMUM LOAD: 220 KG

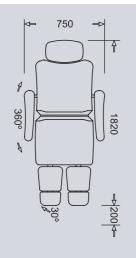


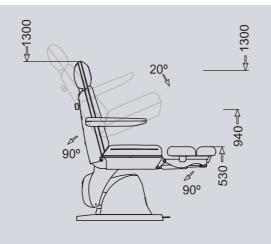

## **DELTA PLUS**


POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 3 A. WEIGHT: 90 KG. MOVEMENT: ELECTROMECHANICAL MOTORS AND ELECTRONICS: LINAK MAXIMUM LOAD: 220 KG

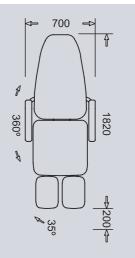


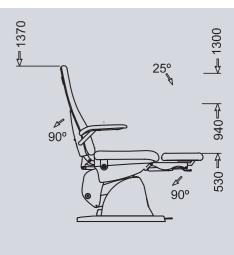

## SIGMA


POWER SUPPLY: 230 V. 50/60 HZ
ABSORBED POWER: 3 A.
WEIGHT: 90 KG.
MOVEMENT: ELECTROMECHANICAL
MOTORS AND ELECTRONICS: LINAK
MAXIMUM LOAD: 220 KG

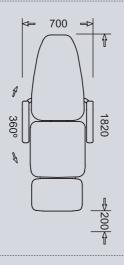


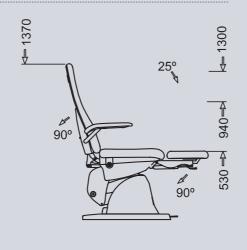

## DELTA


POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 3 A. WEIGHT: 90 KG. MOVEMENT: ELECTROMECHANICAL MOTORS AND ELECTRONICS: LINAK MAXIMUM LOAD: 220 KG

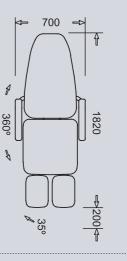


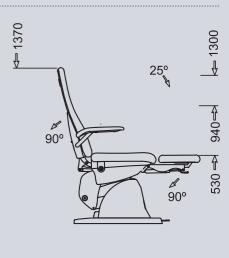

## **OMEGA**


POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 3 A. WEIGHT: 90 KG. MOVEMENT: ELECTROMECHANICAL MOTORS AND ELECTRONICS: LINAK MAXIMUM LOAD: 220 KG

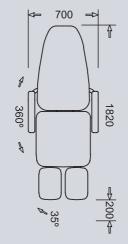


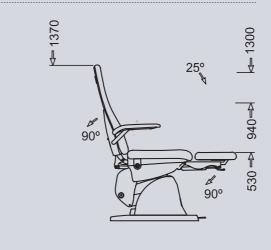

## BETA


POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 3 A. WEIGHT: 90 KG. MOVEMENT: ELECTROMECHANICAL MOTORS AND ELECTRONICS: LINAK MAXIMUM LOAD: 220 KG

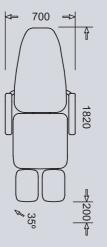


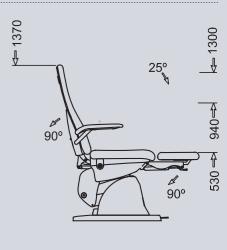

## **ELYSE**


POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 3 A. WEIGHT: 90 KG. MOVEMENT: ELECTROMECHANICAL MOTORS AND ELECTRONICS: LINAK MAXIMUM LOAD: 220 KG

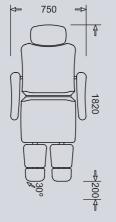


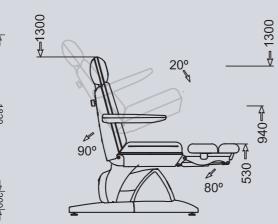

## **PENTA**

POWER SUPPLY: 230 V. 50/60 HZ ABSORBED POWER: 3 A. WEIGHT: 90 KG. MOVEMENT: ELECTROMECHANICAL MOTORS AND ELECTRONICS: LINAK MAXIMUM LOAD: 220 KG

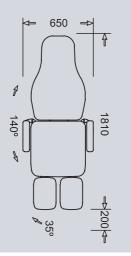

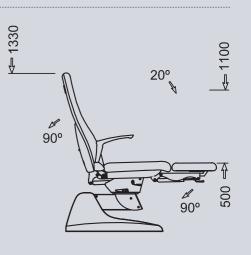





A wide distribution network, with delegations and after-sales service, allows us to reach the entire national territory. Ensuring the best advice of our products and guaranteeing the best service.



#### NAMROL BARCELONA

C/. Progrés, 24-26 08850 Gavà (Barcelona) e-mail: info@namrol.com

### NAMROL VALENCIA

Av. dels Tamarindes, 45 46015 Valencia e-mail: ventasl@namrol.com

### NAMROL MADRID

C/. Puenteáreas, 13 28002 Madrid e-mail: ventasm@namrol.com

### NAMROL SEVILLA

C/. Juan Sebastián Elcano, 15 41011 Sevilla e-mail: ventasa@namrol.com

### NAMROL OVIEDO

Avda. Pedro Masaveu, 9 33007 Oviedo e-mail: comercial@namrol.com



